


THIS DRAWING IS PRELIMINARY AND NOT TO BE USED FOR CONSTRUCTION PURPOSES



EQUIPMENT LIST

ITEM	QTY	DESCRIPTION
01	1	Harmony LED with: 24" HD Screen, LED 785 Lighthead, LED 585 Lighthead, Hub Spacer for Future Lighthead, Ambient Light Accessory and Harmony HD Camera on Auxiliary Mount Arm
02	1	Harmony LED Lighting System Control Center – Standard
03	1	Harmony Remote Power Module (may be mounted in ceiling)
04	1	Vein Harvesting: Tandem Mount with Harmony EMS DF80 with: 40" Supply Column and Harmony LED with: 24" HD Screen, LED 585 Lighthead, and Hub Spacer for Future Lighthead
05	1	Perfusion: Tandem Mount with Harmony EMS DA60 with Supply Head and Harmony LED with: 24" HD Screen and Hub Spacer for Future Lighthead
06	1	Anesthesia and Surgical Equipment: Tandem with Harmony EMS DF80 with 40" Supply Column and DA70 with Supply Head
07	1	STERIS Cmax Surgical Table
08	4	STERIS 42" HD LCD Monitor – Wall Mounted
09	1	STERIS Nurse Documentation Station with: Harmony Integrated OR
10	1	Wall mounted HD Room Camera with pan/tilt/zoom capability

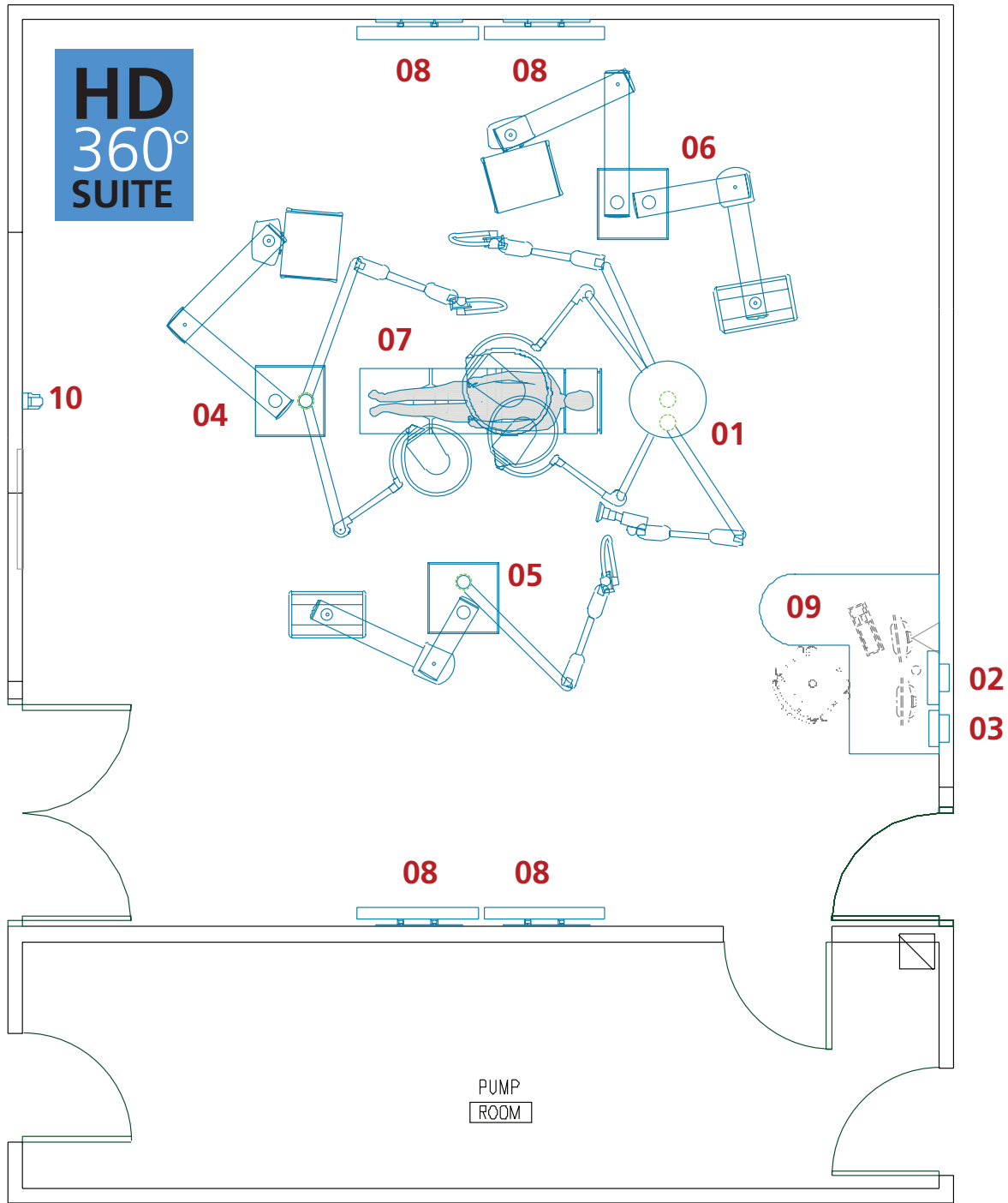

Purpose of the room:

Open Cardiovascular and Thoracic Surgery requires significant lighting, considerable support equipment, the ability to shift to closed procedures, and the pertinent imaging and information at the point of care supported by OR Integration, HD cameras, HD monitors and ceiling support of as much equipment as possible.

Design criteria:

Illumination: Two lightheads for chest (01), one lighthead for vein harvesting (04) and use at the chest

Monitors: Three HD Widescreen Flat Panel Monitors for viewing of pertinent data and archived and real-time images along with four 42" HD Wall Monitors for vitals and other imaging

Perfusion equipment: Perfusion services delivered by dedicated Equipment Management System (EMS) (05)

Endoscopy equipment: EMS (06 or 04) for endoscopic equipment

Anesthesia and surgical equipment: Dual-purpose EMS (06) for anesthesia and for surgical equipment such as defibrillator, headlamp, Chest Bovie, etc.

Design flexibility:

Left or right leg vein harvest: flexible positioning of the EMS unit, lights, and monitor to accommodate either set-up. In the event that perfusion will be located patient RIGHT, you could configure additional utilities on the endoscopic EMS (04) or on the surgeon support EMS (06).

Closed (MIS) cardiovascular procedures: Dual-purpose EMS (06) has potential to hold endoscopic setup to use this room for either open or closed procedures.



CARDIO OPEN – SOME MIS
644 SQ. FT.

HD360 SUITE
EQUIPMENT LAYOUT



STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 • USA
440-354-2600

DR mkf/tmm

CK

DATE 10 April 08

SCALE 1/4" = 1'-0'

H- 9593

B- CVOR1



HD360° SUPPORT

**HD
360°
SUITE**

HD360° Support starts when your room is merely a dream and continues years after your beautiful new space is up and running. We have a dedicated team of Project Design Managers and Systems Planners with over 40 years of experience in designing and planning healthcare facility spaces, and we will help you create the HD360° Suite that meets your needs now and into the future.

We also have the largest field service network in the industry (over 780 strong), who are committed to complete customer satisfaction. Our main priority is to have every STERIS service exceed your expectations.

To see clearly what HD360° total room solutions can do for you, contact your STERIS sales professional or visit www.steris.com/HD360.

STERIS OFFICES WORLDWIDE

Benelux	32 2 523 2488	Japan	81 3 5521 2104
Canada	800 661 3937	Korea	82 2 517 1517
China	86 21 6137 1166	Latin America	800 884 9550
France	33 0 2 38 70 83 50	Malaysia	6 03 7954 9822
Germany	49 2203 890 6969	Singapore	65 68 41 7677
Greece	30 210 6800848	Spain	34 91 658 5920
India	91 33 2367 5150	Switzerland	41 32 376 0200
Italy	39 22 130341	United Kingdom	44 1256 840400

Document # M3217EN.2008-06, Rev. B
Childers Printed 07/2008, 7500

©2008 STERIS Corporation.
All rights reserved. Printed in USA.



STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 • USA
440-354-2600 • 800-548-4873
www.steris.com